

# NIR Spectrometer Ultra-High Optical Resolution

0.04nm, 790 – 880 nm, USB



DATASHEET

[Return to the Webpage](#)



## Features

- Volume Phase Holographic Grating
- Collimating Lens Design
- Focusing Lens Design
- Diffraction System Optical Design
- Wavelength Calibration Algorithms
- Precise Spectrometer Calibrating Techniques
- High Sensitivity
- USB

## Applications

- Laboratory use
- Testing
- Instrumentation
- Optical Coherence Tomography (OCT)

The NRRS Series Spectrometer, based on TE-cooled Silicon detectors coupled with a high throughput transmission grating, is optimized for NIR spectroscopic measurements requiring an exceptional high resolution, signal-to-noise ratio, and high dynamic range across the 790 - 880  $\mu\text{m}$  spectral range. This system leverages innovative, patent-pending scanning technology, offering significant advantages: 1) Unmatched low cost; 2) Industry-leading sensitivity with deep cooling to  $-40^{\circ}\text{C}$ ; 3) Extended spectral coverage beyond traditional spectrometers; 4) Low power consumption; 5) Integrated MEMS chopper; 6) High-resolution performance. The NRRS series comes standard with a USB interface, power supply, and software support includes SDK examples, DLLs for custom application development, and Windows-based spectral acquisition and analysis tools. The NRRS is well suited for RAMAN applications.

## Specifications

Parameter	Min	Typical	Max	Unit
Wavelength Range	790		880	nm
Camera Resolution		4000		pixels
Bit Depth		10, 11, 12		bits
Response non-linearity		$\pm 1$		%
Photo response non-uniformity		0.5		%
Integration Dead Time (in maximum exposure time)	0.01	0.6	2	s
Quantum Efficiency (QE) @850nm		54 %		
Analog gain / Digital gain		x1 x2 and x4 / x1 to x7.996		
Dynamic Range		69		dB
SNR			51	dB
Optical System Characteristics		f/#: 3.6, NA: 0.14 Focal length(R1-R2): 60 -89 @840nm		
Optical Design		T-T-T fully transmissive Czerny-Turner light path		
Grating		1800 lp/mm VPH @840nm		
Input Slit Type		Recommended - 5 $\mu\text{m}$ single mode fiber, FC/PC		
Fiber Optic Interface		FC/PC		
Pixel Resolution	0.035		0.04	nm
Optical Resolution	0.03	0.04	0.07	nm
Operating Temperature	5		40	$^{\circ}\text{C}$
Storage Temperature	-30		70	$^{\circ}\text{C}$
Relative Humidity		0% - 85% non-condensing		
Data transfer interface		USB 3.0		

**Note:** The specifications provided are for general applications with a cost-effective approach. If you need to narrow or expand the tolerance, coverage, limit, or qualifications, please [click this link](#):

**Legal notices:** All product information is believed to be accurate and is subject to change without notice. Information contained herein shall legally bind Agiltron only if it is specifically incorporated into the terms and conditions of a sales agreement. Some specific combinations of options may not be available. The user assumes all risks and liability whatsoever in connection with the use of a product or its application.

Rev 05/07/25

© Photonwares Corporation

**P** +1 781-935-1200

**E** [sales@photonwares.com](mailto:sales@photonwares.com)

**W** [www.agiltron.com](http://www.agiltron.com)

Information contained herein is deemed to be reliable and accurate as of the issue date. Photonwares reserves the right to change the design or specifications at any time without notice. Agiltron is a registered trademark of Photonwares Corporation in the U.S. and other countries.


# NIR Spectrometer Ultra-High Optical Resolution

0.04nm, 790 – 880 nm, USB



## DATASHEET

### Mechanical Dimensions (mm)

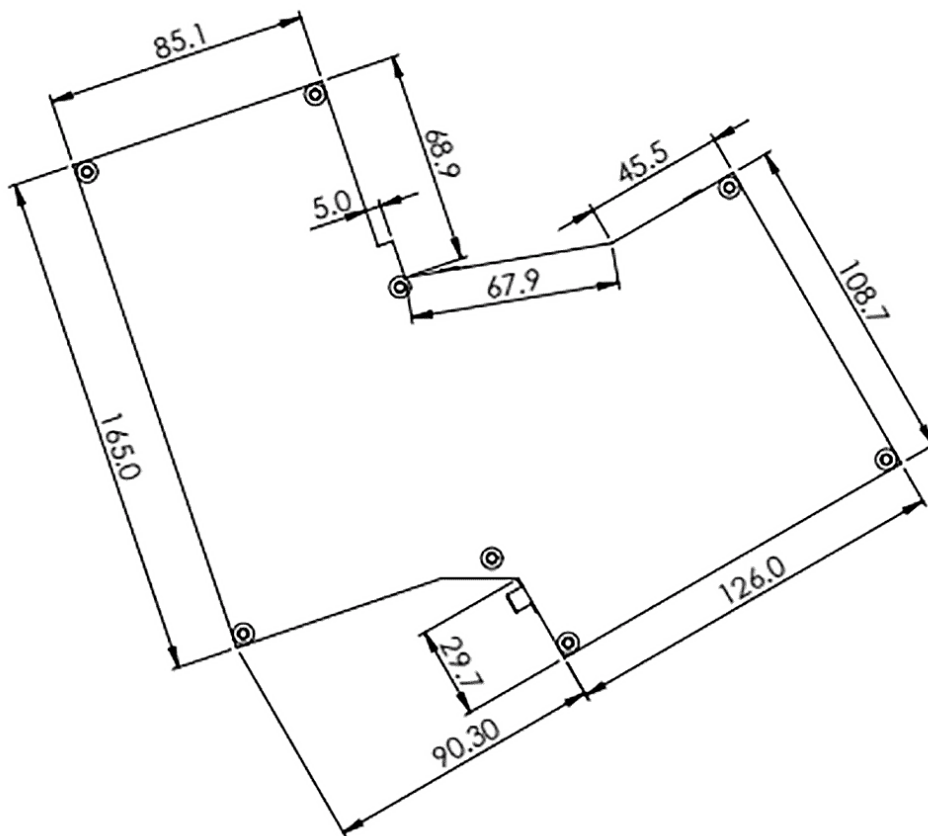

\*Product dimensions may change without notice. This is sometimes required for non-standard specifications.

# NIR Spectrometer Ultra-High Optical Resolution

0.04nm, 790 – 880 nm, USB



## DATASHEET

### Ordering Information

	1	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>	<input type="checkbox"/>
Prefix	Type	Wavelength	Input Optical Power	Cooling *	Resolution **	Shutter	Chopper	Connector
HRRS-	Raman = 1	790-880nm = 1 Special = 0	Standard = 1 High Power = 2	Non = 1 -5°C = 2 -20°C = 3 -40°C = 5	1nm = A1 0.5nm = 50 0.05nm = 05 0.04nm = 04 Special = 00	Non = 1 Yes = 2	Non = 1 Yes = 2	SMA905 = 1 FC/PC = 2 SC/PC = 4 ST/PC = 6 Special = 0

\* Non cooling is low cost for strong light measurements. At -5°C:Noise is reduced by about 4×, improving performance in low-light applications. At -40°C:Noise is reduced by about 16×, enabling high-sensitivity measurements, such as weak signal detection in spectroscopy or astronomy.

\*\* Low resolution high sensitivity.